

LIGHT CABLE TESTER

NORUCO

- EASY TO USE

- VERY INTUITIVE INTERFACE

- EXACT MEASUREMENTS

- RESULTS IN ONLY 5 SECONDS

- HIGH QUALITY COMPONENTS

- REPETITION OF THE TEST AS MANY TIMES AS IT WANTS


The **Noruco NC10** tester gives the possibility to check the status of the light cables. For this, it uses high quality components, giving a very accurate and safe measurement result.

You can use the tester with any type of connection from both the endoscope and the cold light source and from any manufacturer, thus avoiding unnecessary substitutions and keeping our light transmission systems always in optimal condition.

The equipment facilitates the work of both the medical / surgical personnel and the maintenance personnel.



NC10

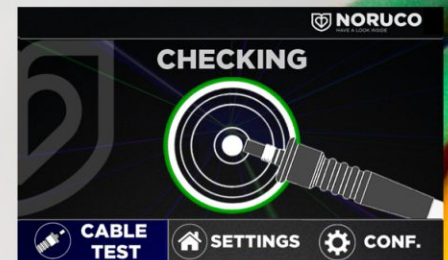

The **Noruco NC10** tester compares the luminous flux between both ends of the cable and processes the difference thanks to an accurate calibration.



We set various levels of light cable status:

100%-80%



Optimum

79%-60%



Valid

59%-0%



Not valid